

FORTGESCHRITTENEN PRAKTIKUM II

---

# Optical pumping

---

11.04.2016

Benjamin Winkelmann  
Peter Spalthoff

Tutor: Dominik Schomas

# Contents

List of Figures	II
<b>1 Goal of the experiment</b>	<b>1</b>
<b>2 Physical principles</b>	<b>1</b>
2.1 Hyperfine structure and Zeeman splitting . . . . .	1
2.2 Optical pumping . . . . .	3
2.3 Relaxation processes . . . . .	4
2.4 Larmor precession of the spin . . . . .	5
<b>3 Experimental set-up</b>	<b>6</b>
<b>4 Characterization of the laser diode</b>	<b>7</b>
<b>5 Spectroscopy of the hyperfine structure</b>	<b>8</b>
5.1 Set-up and procedure . . . . .	8
5.2 Data analysis . . . . .	9
<b>6 References</b>	<b>12</b>

## List of Figures

2.1	Hyperfine structure of Rubidium . . . . .	2
2.2	Hyperfine structure energies . . . . .	3
2.3	Optical pumping . . . . .	4
3.1	Basic experimental set-up . . . . .	6
4.1	Characteristic curve of the laser diode . . . . .	7
5.1	Experimental set-up for the Calibration . . . . .	8
5.2	Bend in modulating current . . . . .	9
5.3	Etalon peaks . . . . .	10
5.4	Linear fit on etalon peak positions . . . . .	11
5.5	The HFS-spectrum . . . . .	13

## 1 Goal of the experiment

In this experiment, the process of optical pumping will be used to precisely measure properties of Rubidium atoms such as the hyperfine constant  $A$  via absorption measurements. In addition to that, relaxation times of the induced pumped states as well as external magnetic fields will also be measured by observing the effect of magnetic fields, applied through Helmholtz coils, high frequency radio waves and variations in the laser intensity.

## 2 Physical principles

### 2.1 Hyperfine structure and Zeeman splitting

This section is based on the detailed elaborations in [1].

The fine structure levels of the atomic spectrum, which splits the basic levels into sub-levels due to spin-orbit interaction, can be shown to be split into even finer levels, whose energetic distances are roughly three orders of magnitude smaller than those of the fine structure. This is called the 'hyperfine structure' and is mainly caused by the interaction of the nuclear magnetic dipole and quadrupole moment and the magnetic field of the shell electrons. Its structure for the two Rubidium isotopes that are used in this experiment can be seen in figure 2.1.

As the nucleus is charged and, expressed as the nuclear spin  $\vec{I}$ , has angular momentum, it also has a magnetic moment, which is  $\vec{\mu}_I = \frac{g_I \mu_K}{\hbar} \vec{I}$ , where  $g_I$  is the g-factor of the nucleus and  $\mu_K$  is the nuclear magneton.

With the total angular momentum of the electrons  $\vec{J}$ , the total angular momentum of the atom can be written as

$$\vec{F} = \vec{J} + \vec{I}, \quad |I - J| \leq F \leq I + J \quad (2.1)$$


The energy difference between hyperfine structure levels can then shown to be

$$\Delta E_{HFS} = -\vec{\mu}_I \cdot \vec{B}_J = \frac{A}{2}(F(F+1) - J(J+1) - I(I+1)) \quad (2.2)$$


where  $A = \frac{g_I \mu_K B_J}{\sqrt{J(J+1)}}$  is the hyperfine constant. Neighboring levels thus have an energy difference of

$$\Delta E_{HFS}(F+1) - \Delta E_{HFS}(F) = A(F+1) \quad (2.3)$$

This structure for the rubidium isotopes used in this experiment can be seen in figure 2.2.

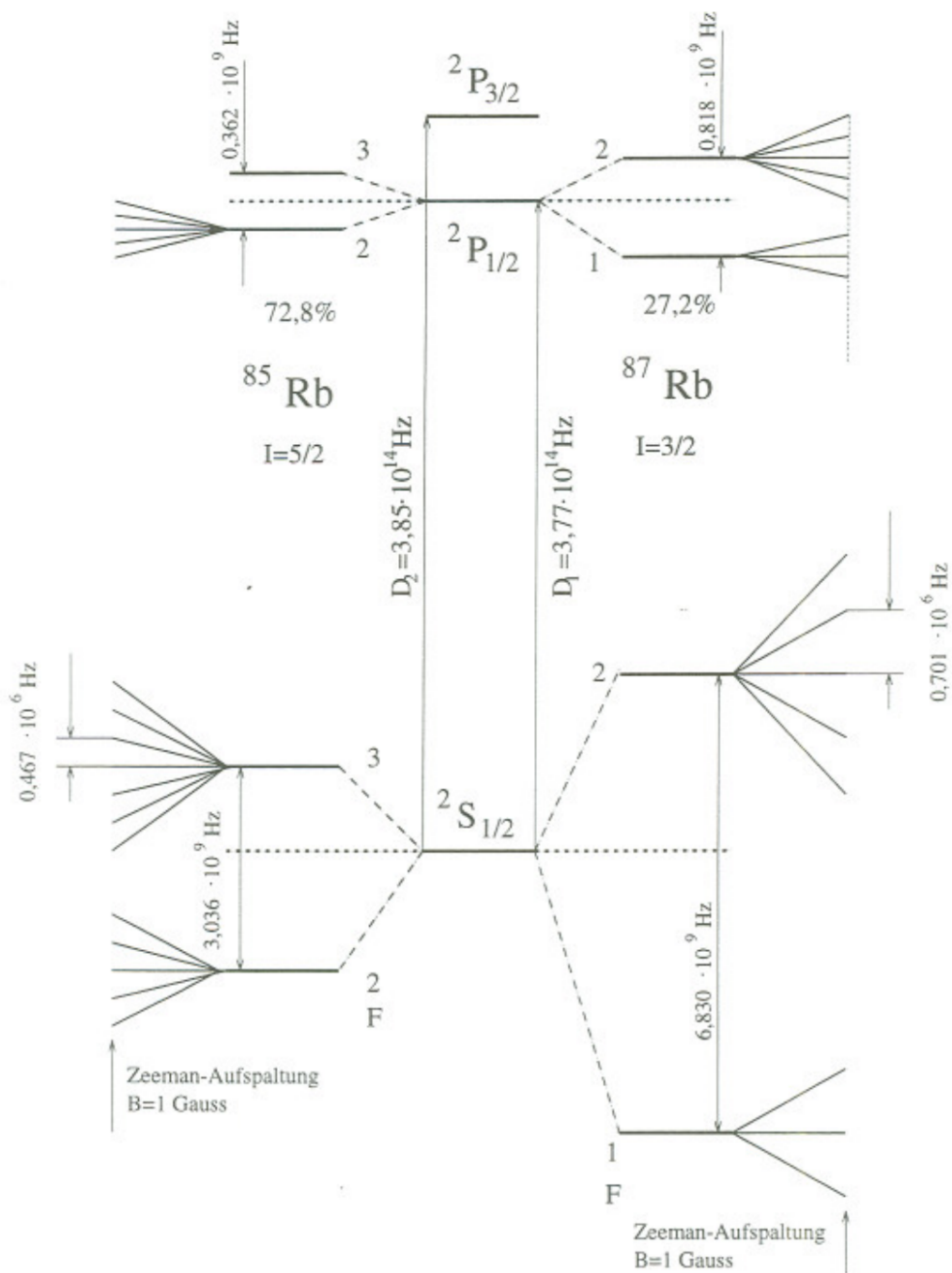

**Figure 2.1:** The hyperfine structure of the two isotopes of Rubidium used in the experiment. The hyperfine levels in turn are split due to the Zeeman effect caused by an external field of  $B = 1 \text{ G}$ . [1]

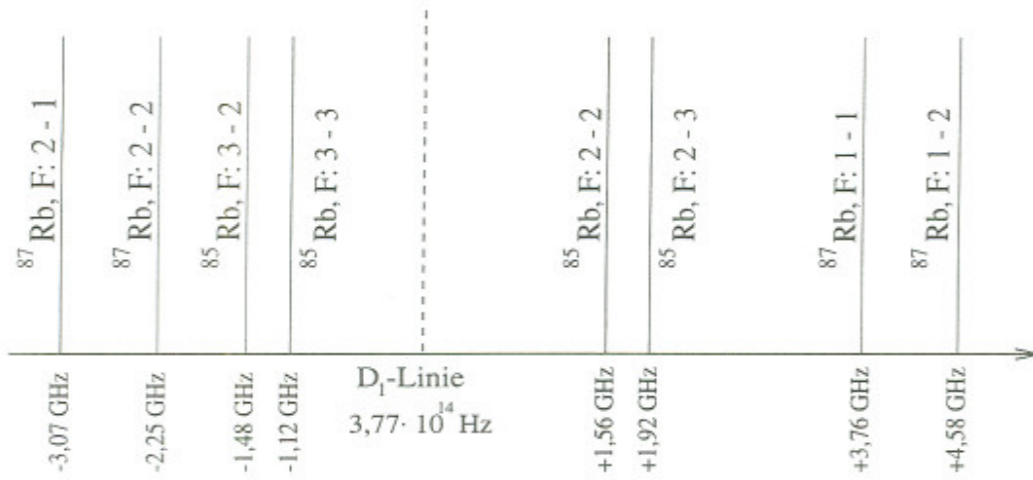

**Figure 2.2:** Energies of a set of hyperfine levels of the two isotopes used in the experiment. [1]

These hyperfine levels can in turn be split into  $2(F + 1)$  sub-levels in the presence of an external magnetic field. The according quantum number is  $-F \leq m_F \leq F$ . As long as the magnetic field is weaker than the spin-orbit coupling or, in other terms,  $g_J \mu_B B_0 \ll A$ , this is called the Zeeman effect. The effect for larger fields, where the spin-orbit coupling is disrupted, is called Paschen-Back effect.

For the Zeeman effect, the energy difference of the levels is

$$\Delta E_{Zeeman} = \frac{g_J}{2(I + \frac{1}{2})} \mu_B B_0 \quad (2.4)$$

## 2.2 Optical pumping

In general, pumping refers to constantly transferring electrons into higher energy levels until significantly more electrons are in the higher than in the lower state. This is called population inversion.

In the case of this experiment, this is done using a laser diode. As the goal is to examine magnetic fields using the Zeeman splitting, a way must be found to create population inversion within a single non-degenerate hyperfine structure level. Normally, electrons are equally distributed between said levels.

The selection rules for transitions

$$\begin{aligned} \Delta F &= 0, \pm 1 & (F = 0 \leftrightarrow F = 0) \\ \Delta m_F &= 0, \pm 1 \end{aligned} \quad (2.5)$$


allow for a convenient way to change that. If only  $\sigma^+$ -polarized light is used, only transitions with  $\Delta m_F = +1$  are caused. Since the following decay is random within the bounds of the transition rules, the laser will pump all electrons into the  $^2S_{1/2}$  state with

$m_F = +2$ ,  $F = 2$  for  $^{87}\text{Rb}$  and  $m_F = +3$ ,  $F = 3$  for  $^{85}\text{Rb}$ . Figure 2.3 illustrates this for two exemplary transitions.

Mathematically, this process can be described as

$$\left(\frac{dn}{dt}\right)_P = \frac{N - n}{T_P} \quad (2.6)$$

where  $n$  is difference of the levels in the two-level system,  $N$  the overall number of atoms in the system and  $T_P$  the characteristic pumping time of the system.

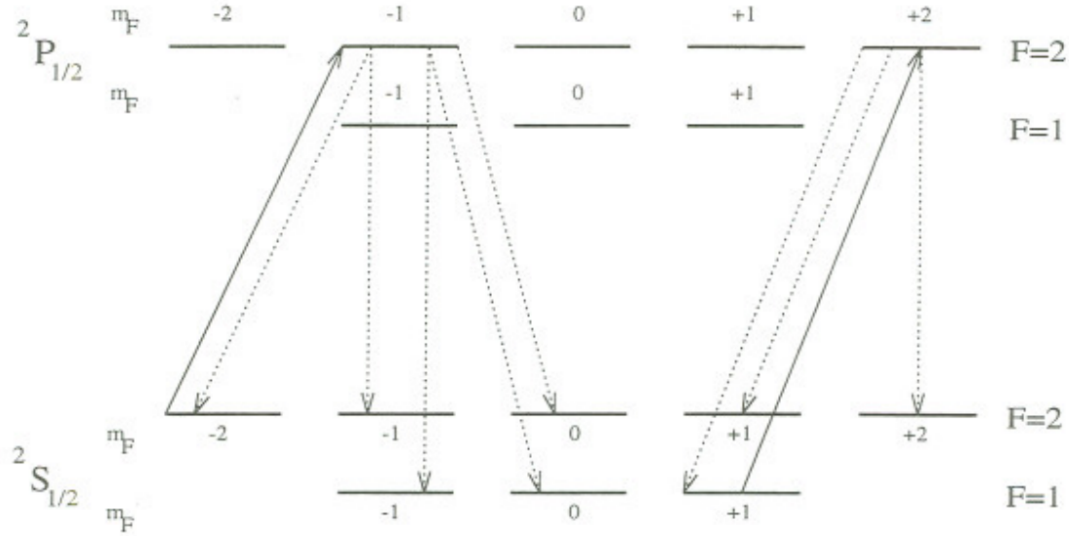

**Figure 2.3:** Optical pumping for  $^{87}\text{Rb}$ . The  $\sigma^+$  polarized light can only cause transitions with  $\Delta m_F = +1$ , thus achieving the desired pumping effect.

### 2.3 Relaxation processes

The desired pumping effect is counteracted mainly by three relaxation effects:

**Diffusion to the wall:** Upon hitting the glass containment, the rubidium atoms may lose their polarization. This process is inhibited by a buffer gas, which limits the mean free path of the atoms.

**Collisions with the buffer gas:** The rubidium atoms may lose their polarization due to collisions with the buffer gas. The cross section of this event depends highly on what kind of gas is used. Best results are achieved with noble gases - in the case of this experiment, krypton was used.

**Spin exchange** When rubidium atoms collide, they may interchange their spins. While the overall polarization is preserved, the decoupling of nuclear and electron spins lead to a faster relaxation time. Much more detailed elaborations can be found in [3].

Overall, the relaxation can be described by the following differential equation

$$\left(\frac{dn}{dt}\right)_R = -\frac{n}{T_R} \quad (2.7)$$

where  $T_R$  is the characteristic relaxation time. A value of  $T_R^{theo} = [6.5]ms$  is given in [1]. The overall process of polarization orientation is thus the sum of equations 2.6 and 2.7:

$$\left(\frac{dn}{dt}\right)_O = \left(\frac{dn}{dt}\right)_P + \left(\frac{dn}{dt}\right)_R = \frac{N}{T_P} - n\left(\frac{1}{T_P} + \frac{1}{T_R}\right) \quad (2.8)$$

The solution of this equation is an exponential

$$n(t) \propto e^{-\frac{t}{\tau}} \quad (2.9)$$

where  $\tau = \frac{1}{T_P} + \frac{1}{T_R}$ .

## 2.4 Larmor precession of the spin

If the ensemble is polarized along a certain magnetic field and one component of said field is suddenly set to zero, the polarization precesses around the remaining field. The precession frequency is


$$f_L = \frac{g_F \mu_B}{h} \cdot B =: \alpha \cdot B \quad (2.10)$$

where  $g_F$  are the Landé factors for the rubidium isotopes, quantified by Baur [1] as  $g_F(^{85}Rb) = 1/3$  and  $g_F(^{87}Rb) = 1/2$ . The proportionality constant between the frequency and the remaining magnetic field thus is

$$\alpha(^{85}Rb) = 4.665 \text{ kHz}/\mu\text{T} \quad \alpha(^{87}Rb) = 6.998 \text{ kHz}/\mu\text{T} \quad (2.11)$$



### 3 Experimental set-up

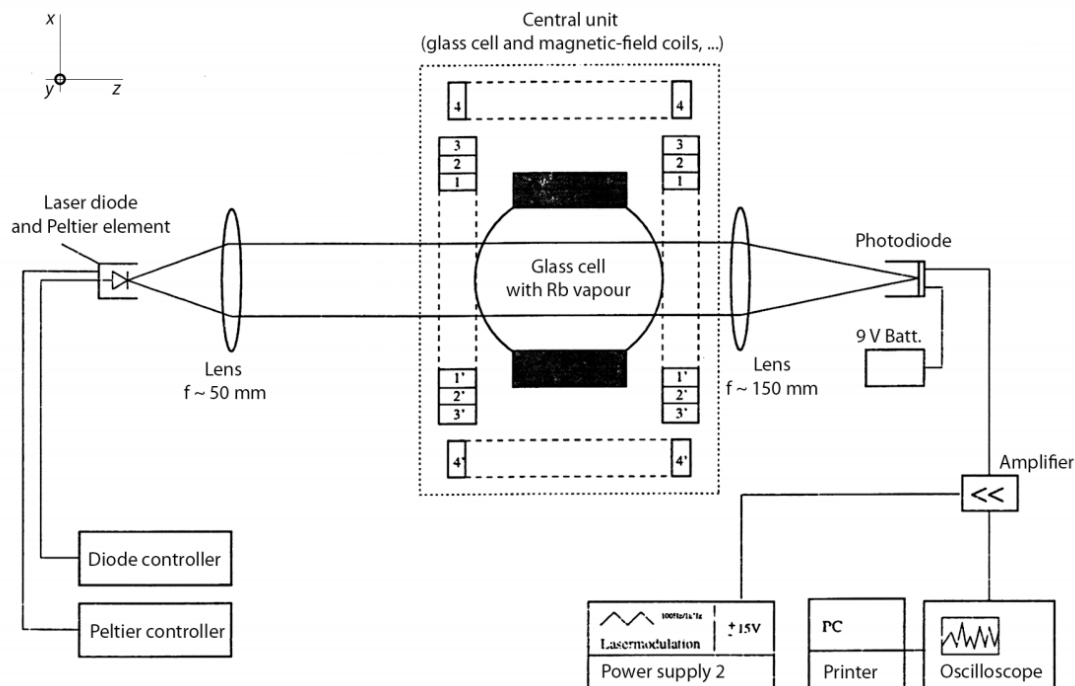

**Figure 3.1:** Basic set-up of the experiment. For the various tasks, parts can be placed onto the optical bench. The glass cell can be removed from the central unit. [2] (modified)

At the core of the experimental set-up is the glass cell with the Rubidium vapor and the buffer gas. It can be taken in and out of the central unit, which in turn houses four sets of Helmholtz coils. Another such pair is directly attached to the casing of the glass cell, along with a radio frequency generator and an appropriate frequency measuring device.

**A laser diode** provides coherent light in the energy range needed to pump the desired hyperfine state of the Rubidium atoms in the glass cell. Laser diodes send out linearly polarized light with a small spectral width. Frequency and intensity vary with the temperature of the diode as well as the current running through it, which is why the diode is kept at constant temperature using a Peltier element. For the diode to start emitting light, a certain current threshold has to be reached. From then on, the intensity depends linearly on the current if temperature is kept constant. However, mode jumps at certain currents disturb the linearity. Mode jumps occur when the number of standing waves in the resonator changes and measurements need to be taken in areas that do not include such jumps.

The beam is collimated by a lens before passing through other optical elements and,

after passing through the central unit, is refocused onto a photo-diode. This can be seen in figure 3.1. The output of the diode is amplified and can then be observed on an oscilloscope, which in turn can be read out by a computer to produce analyzable data.

The set-up varies greatly from one part of the experiment to another and will thus be explained in detail in the appropriate sections.


## 4 Characterization of the laser diode

For later measurements, it is important to determine the range of supply current in which the diodes intensity increases linearly without mode jumps occurring. The gas cell is taken out of the central unit for this part of the experiment.

After turning on the peltier element, a few minutes should pass before measurements are started to allow the diode to thermalize. Measurements were taken at  $T = 34.3^\circ$ .

Since the photo diode saturated for laser diode currents upwards of  $I_L = 65$  mA, a neutral filter (D2,6, see figure ??) is used to limit intensity so that the photo diode barely does not reach saturation.


The intensity of the diode is now measured at supply currents between 0 – 90 mA. The results can be seen in figure 4.1. The threshold current is roughly 51.6 mA, followed by the linear domain until mode jump occur at around 72 mA to 82 mA.

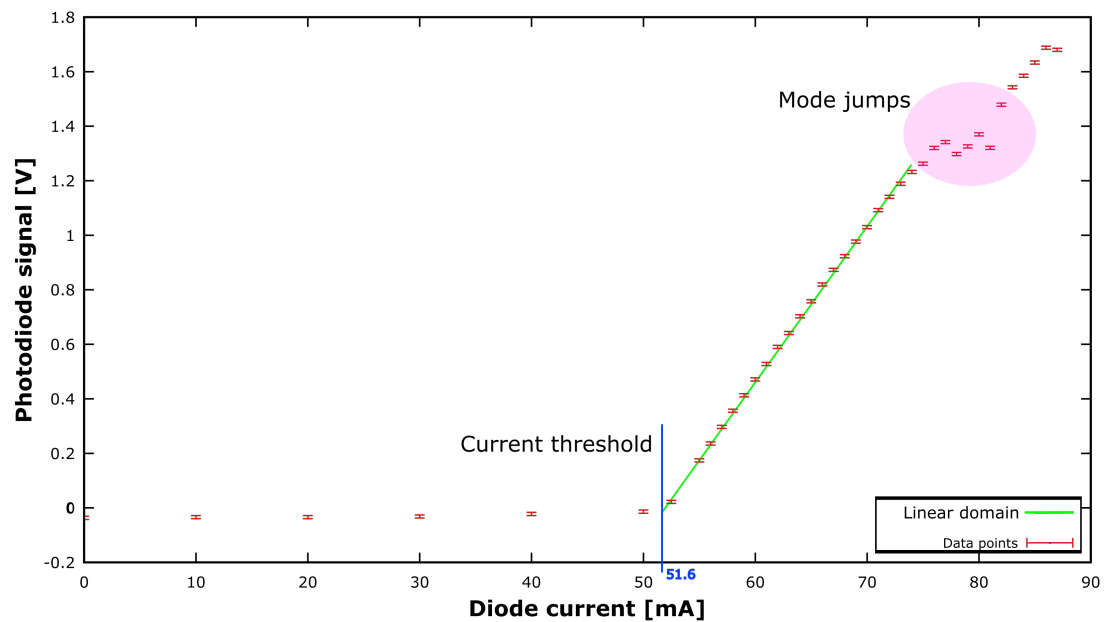

**Figure 4.1:** The characteristic curve of the laser diode for supply currents up to 90 mA. A mode jump can clearly be seen in the upper right corner.

## 5 Spectroscopy of the hyperfine structure

### 5.1 Set-up and procedure

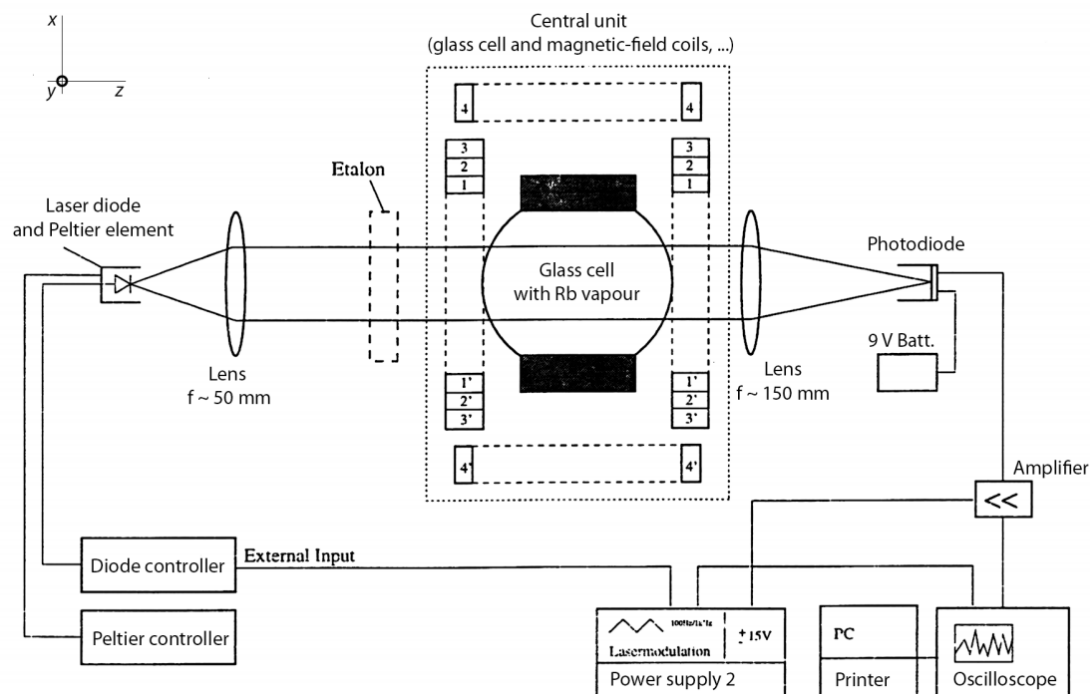


**Figure 5.1:** Experimental set-up for the calibration measurements. An etalon is placed after the left lens and can be removed for the HFS-spectrum measurement.[2]

#### Laser scan-rate

A Fabry-Pérot-Interferometer (etalon), which is a set of parallel, almost completely reflecting surfaces facing each other, is used to gauge the diode. The intensity at the photo diode will be drastically dampened unless the laser wavelength allows for constructive interference to occur between the surfaces. The wavelengths at which light passes the interferometer are thus equidistant.

The diode current is then modified with a saw-tooth voltage, which is also used to trigger the oscilloscope (see figure 5.1). It is important to measure the scan-rate in a range of diode currents where the intensity response is linear and where the HFS-spectrum is actually visible (see next section). To get a solid calibration, at the very least 3 etalon peaks are necessary. These two factors have to be taken into consideration when choosing the measurement interval.

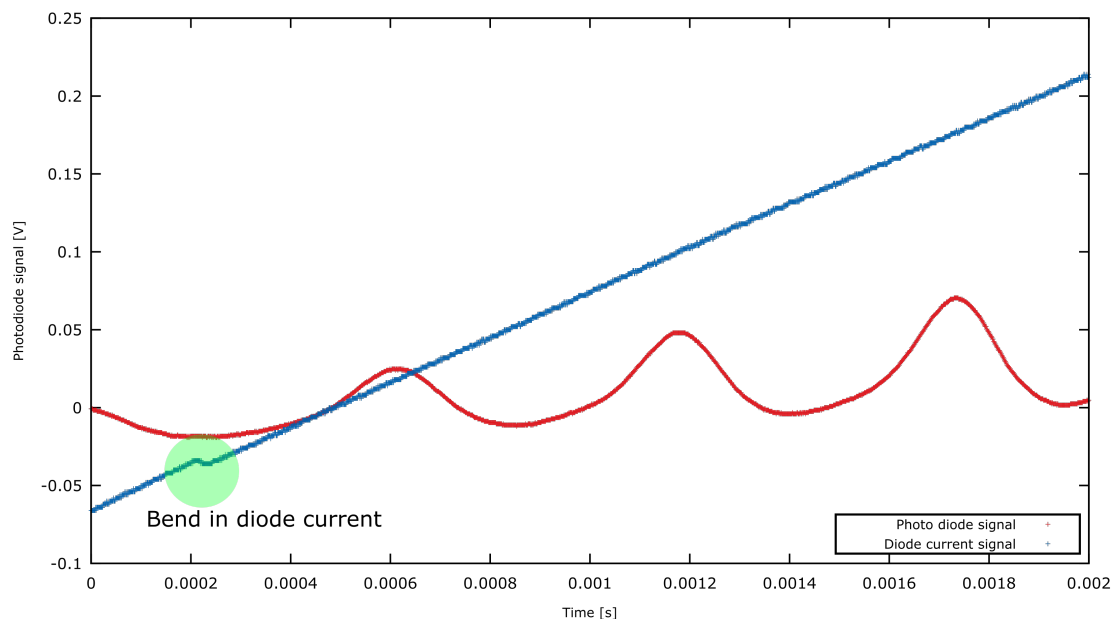

**Figure 5.2:** In the bottom left corner, a clear bent in the diode current can be seen. These measurements were not used, a better current range was chosen instead.

Figure 5.2 shows a bend in the diode current. Such areas need to be avoided when taking measurements, as they would lead to non-equidistant peaks in the etalon spectrum and unnecessarily complicate data analysis. The base current for the measurements was  $I_L = (58.5 \pm 0.1)$  mA, where the uncertainty is the digit error of the multimeter ELTELEC, for which no data sheet was available. The temperature remained at  $T = 34.3^\circ$ . The resulting intensity distribution can be seen in figure 5.3.


### The HFS-spectrum

The etalon is removed for this measurement and the exact same current modulation of the laser current as has been used for the scan-rate measurements will be used here as well. As the intensity increases with increasing diode current, the spectrum would have a linear offset. This can be counteracted by inverting the photo diode signal and adding the diode current signal to it. An additional cosmetic advantage of this is that the peaks are then positive.

## 5.2 Data analysis

### Laser scan rate

In figure 5.3, it has to be noted that the saw-tooth voltage reaches its maximum between the rightmost peak and the one on the left of it, meaning that these two actually belong

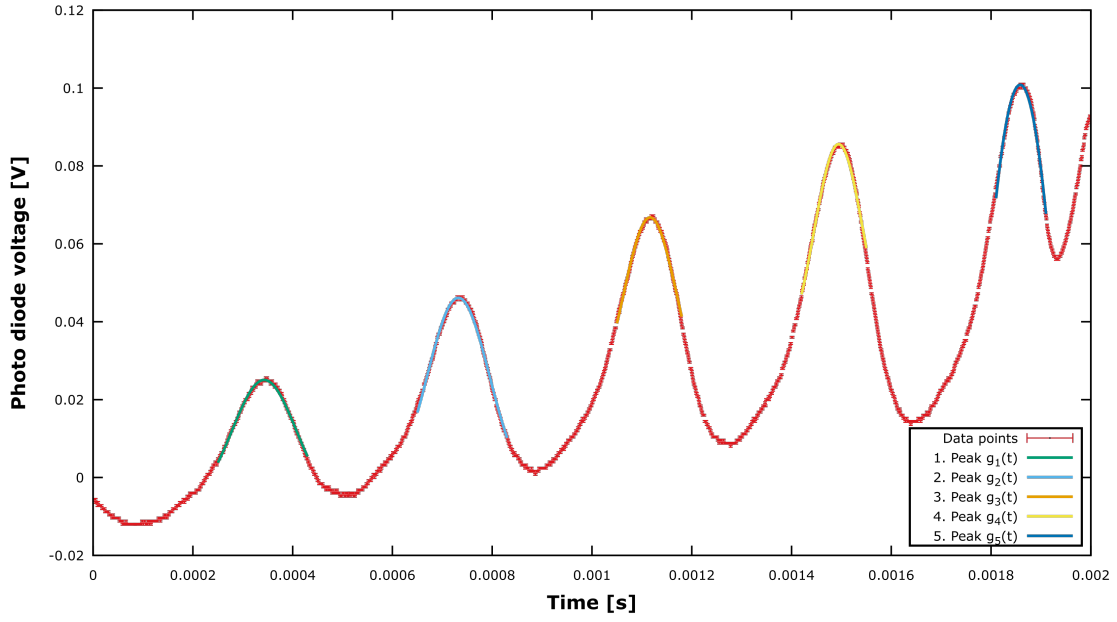

**Figure 5.3:** Results of the etalon calibration measurements. The rightmost peak is actually the same as the one on the left of it, as the saw-tooth reaches its peak between the two.

to the same frequency. The rightmost peak is thus ignored for the calibration. Gauss functions with a constant offset of the form


$$g_i(t) = a_i * e^{-\frac{(t-t_{0,i})^2}{2*\sigma_i^2}} + c_i \quad (5.1)$$

and a linear offset function  $l(t) = b * t + c$  were used. The results for the peak positions  $t_{0,i}$  can be seen in figure 5.4. As the linear fit shows, these peaks are equidistant as expected with the time distance between peaks being  $b = (0.379 \pm 0.003)$  ms. The nominal distance between two peaks is  $\Delta_{FSR} = (9924 \pm 30)$  MHz. With these two values, the scan rate  $r$  of the laser can be calculated as

$$r = (26.17 \pm 0.22) \frac{\text{GHz}}{\text{ms}} \quad (5.2)$$

### HFS-spectrum

Only 6 of the expected 8 peaks (see figure 2.2) are visibly separated here. It is reasonable to assume that the peaks that were not separable are those of the transitions  $^{85}\text{Rb}$  F:3-2 and  $^{85}\text{Rb}$  F:3-3 as well as  $^{85}\text{Rb}$  F:2-2 and  $^{85}\text{Rb}$  F:2-3. This also means that the hyperfine constant for the  $^2S_{1/2}$  of  $^{85}\text{Rb}$  cannot be calculated with these measurements. Gauss functions like those in equation 5.1 were used for the individual peaks. However, for the peaks that overlap, as sum of the two functions was used. For all fits, both a linear

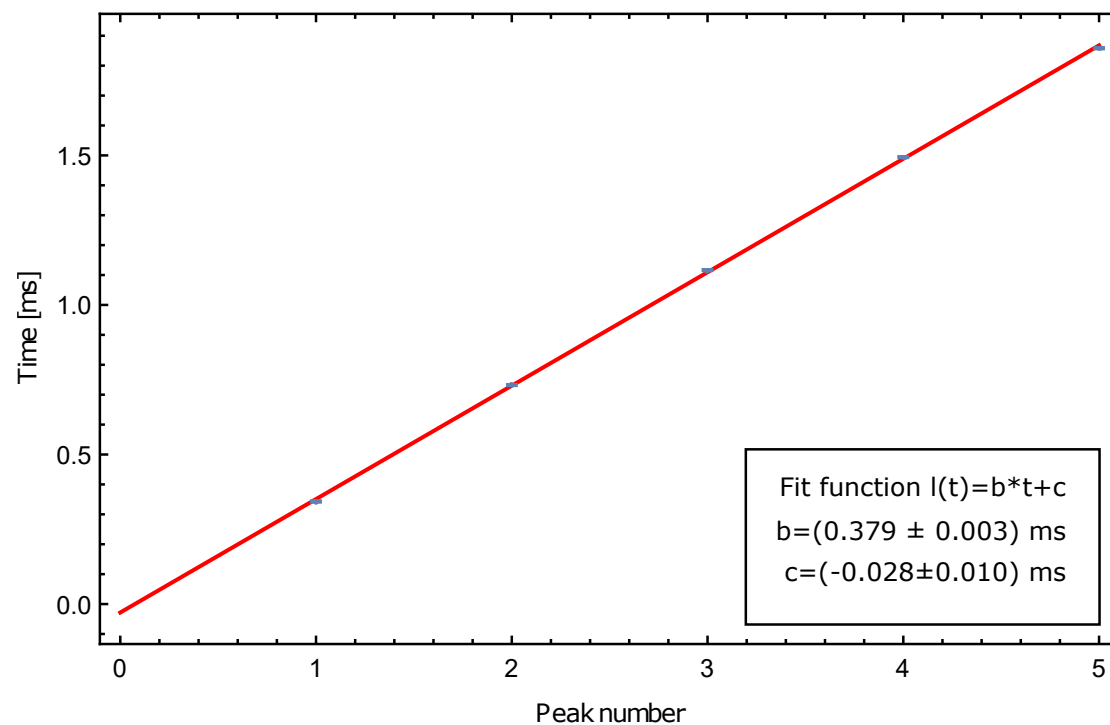

**Figure 5.4:** Linear fit on the etalon peak positions.

function and a constant were added and the total sum fitted to the data. The fit function for the last three peaks for example was

$$I_{4-6}(t) = g_4(t) + g_5(t) + g_6(t) + b \cdot t + c \quad (5.3)$$


Table 5.1 lists the result of the fits as well as the frequency distances to the reference peak, which was arbitrarily chosen to be the first peak.

Peak No	$t_0$ [ms]	$s_{t_0}$ [ $\mu$ s]	$\tau^+\tau^-$	Transition
1	0.104	0.18		$^{87}\text{Rb}$ F:1-2
2	0.129	0.92		$^{87}\text{Rb}$ F:1-1
3	0.215	0.05		$^{85}\text{Rb}$ F:2-2/3
4	0.340	0.03		$^{85}\text{Rb}$ F:3-2/3
5	0.371	0.06		$^{87}\text{Rb}$ F:3-2
6	0.399	0.06		

**Table 5.1:** Fit results of the Gauss fits on the HFS spectrum. The frequency calibration was used to calculate the differences to the reference peak.

## 6 References

- [1] Baur, Clemens. "Einrichtung des Versuchs *Optisches Pumpen mit Laserdioden*." Zulassungsarbeit. Freiburg 1997
- [2] "Instructions for the experiment *Optical Pumping*." Albert-Ludwig University Freiburg. Freiburg 02.2016
- [3] Happer, William. "Optical pumping." Reviews of Modern Physics 44.2 (Apr. 1972): 170-238

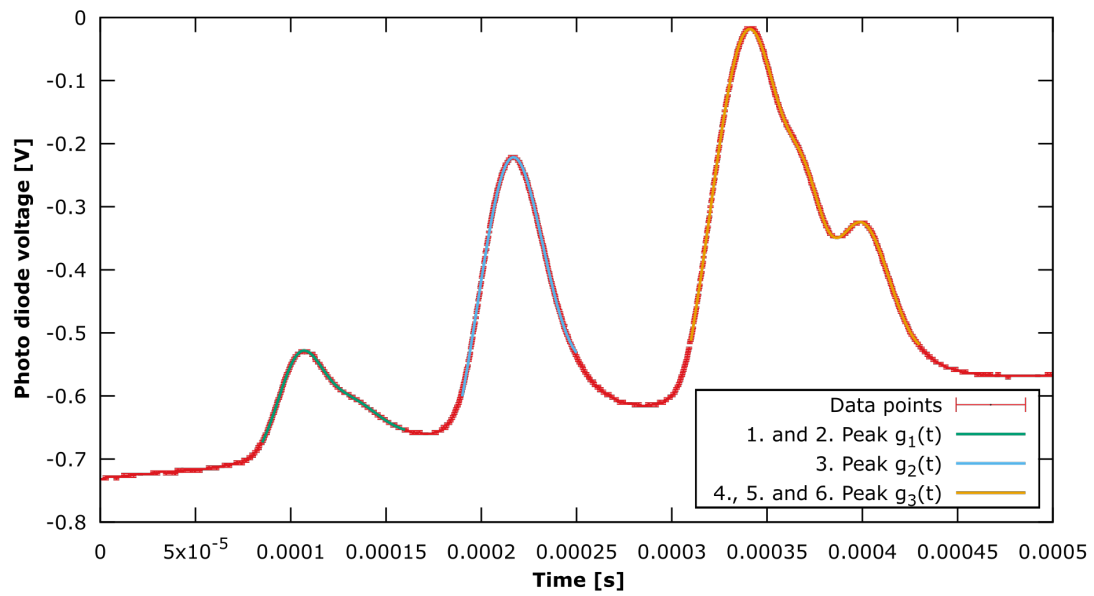

**Figure 5.5:** The HFS-spectrum with fits to the individual peaks. The signal is negative since it was inverted.